

TSIV SA Fixture 03.2010

- Submerged type implant with an internal hex 11° tapered connection structure
- Optimal thread design for realization of optimal SA surface
- Dedicated fixture for use in maxillary sinus and soft bone
- Excellent initial stability in soft bone with smaller threads in the upper section
- Superior self-threading effect with corkscrew thread
- Sharp apex design allowing placement even after D4 bone Ø2.0/3.0mm drilling

Ultra-wide


- Ideal for placement in a fresh extraction socket in the posterior area or for replacing a failed implant
- Optimized apex design for excellent initial stability in a fresh extraction socket or in 3mm from the bottom
- Recommended placement torque : ≤40Ncm
- ※ Fixtures with a diameter of 4.5mm or greater are recommended for the posterior region
- ※ Reducing the speed to 15rpm or lower recommended for placement as the placement speed is too fast for TSIV fixtures due to large thread pitch

NoMount Fixture order code

: fixture product code (ex : TS4S4010S)


Pre-Mounted Fixture (fixture + mount + cover screw) order code

: B + fixture product code (ex : B**TS4S4010S**)

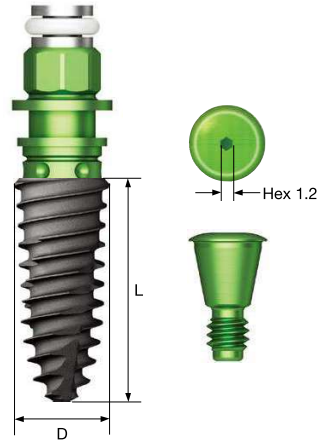


Ultra-wide

D Ø6.0
Hex 2.5

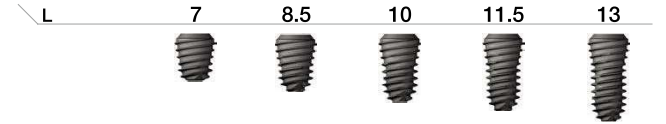


TS4S6007S TS4S6008S TS4S6010S TS4S6011S TS4S6013S

D Ø7.0
Hex 2.5

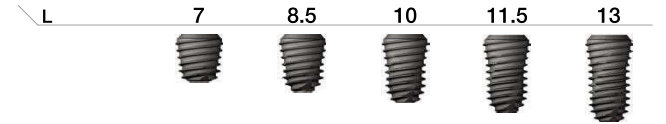


TS4S7007S TS4S7008S TS4S7010S TS4S7011S TS4S7013S

D Ø4.0 Pitch 0.8
Hex 2.5

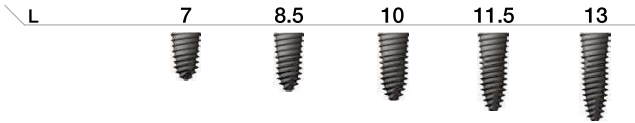


TS4S4007S TS4S4008S TS4S4010S TS4S4011S TS4S4013S

D Ø4.5 Pitch 1.0
Hex 2.5

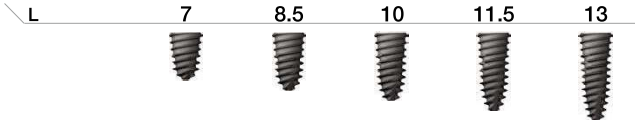

TS4S4507S TS4S4508S TS4S4510S TS4S4511S TS4S4513S

D Ø5.0 Pitch 1.2
Hex 2.5

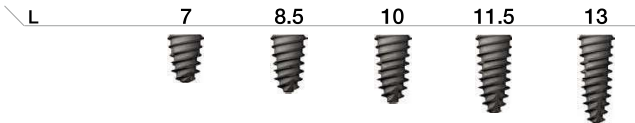

TS4S5007S TS4S5008S TS4S5010S TS4S5011S TS4S5013S

Nominal diameter may differ from actual diameter.